


GEMÜ B600 02-01.H

Tank valve

EN

Assembly instructions

Anbau auf Behälterdeckel

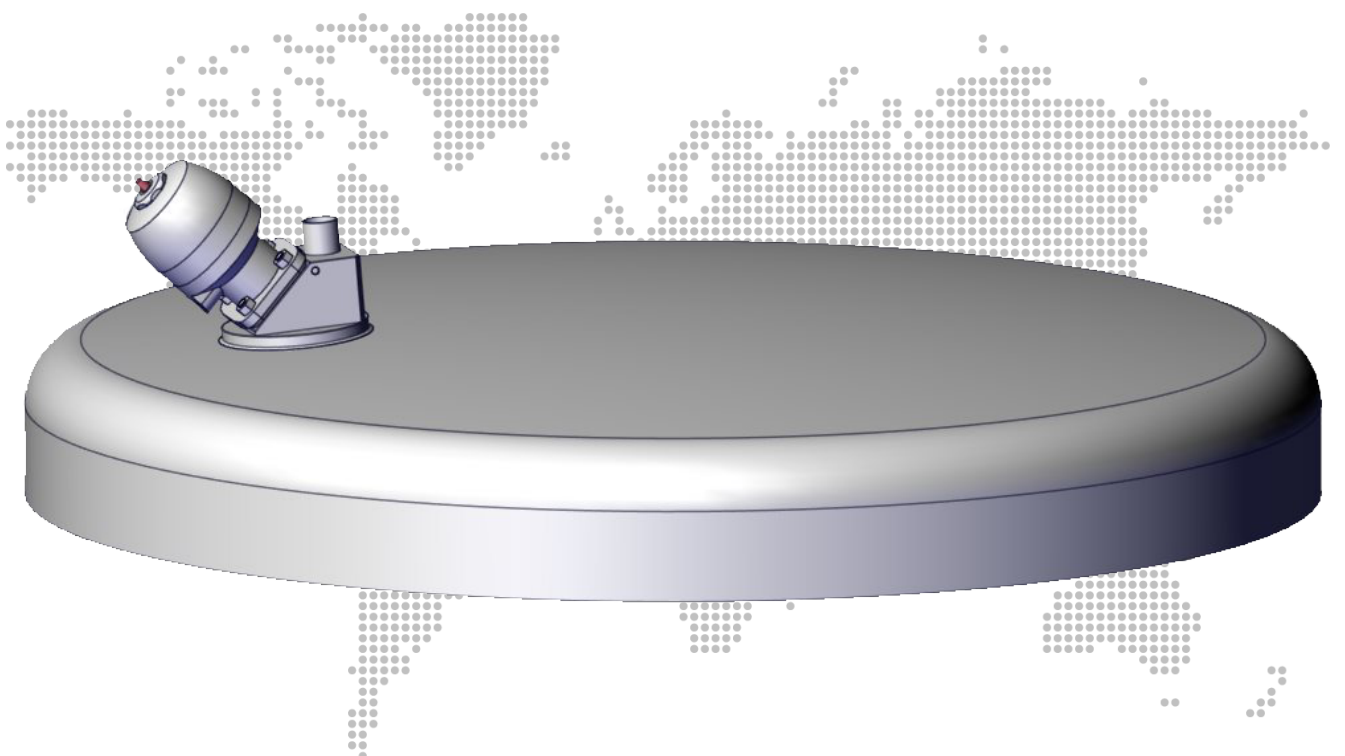


1 Attaching to the tank cover

NOTICE


Valve installed incorrectly

- ▶ Valve does not run completely empty.
- ▶ Residues can cause the medium to become contaminated.
- During installation of the valve body, ensure the inclination angle is correct.

1. Remove the actuator from the valve body.
2. Position the valve body on the tank cover so that the medium flows through the inlet nozzle **1** to the valve seat **2**.

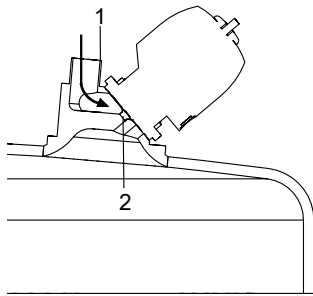

Correct inclination angle

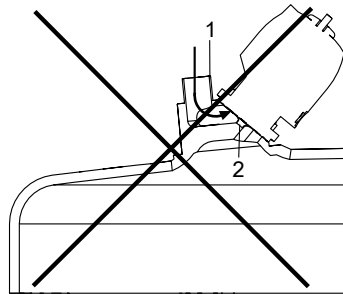

Incorrect inclination angle

3. Weld in the valve body.
4. Fit the actuator to the valve body.



Subject to alteration

02.2019 | 88616768

